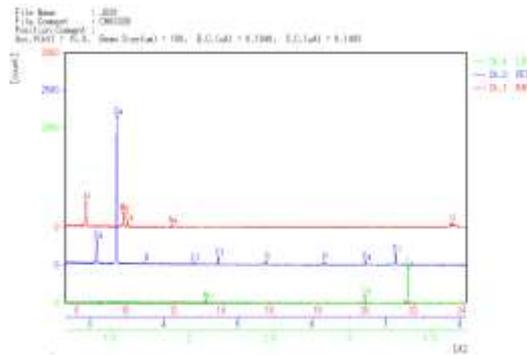


| 茶渋の元素分析 <<株式会社コーエキ>> | | 1009 EPMA/元素分析 | | | | | | | | | | | | | | | | | | | | | | | | |
|----------------------|---|---|----|---|------|----|------|----|-----|----|-----|----|-----|----|-----|---|-----|----|-----|---|-----|----|-----|---|-----|---|
| 試料情報 | 茶碗に付着した茶渋 |  | | | | | | | | | | | | | | | | | | | | | | | | |
| 装置名 | 電子線マイクロアナライザー (EPMA-1610) | | | | | | | | | | | | | | | | | | | | | | | | | |
| 測定条件 | 採取、カーボン蒸着、元素範囲： ₈ O (酸素) ~ ₉₂ U (ウラン) | | | | | | | | | | | | | | | | | | | | | | | | | |
| 備考 | | | | | | | | | | | | | | | | | | | | | | | | | | |
| 測定結果 | <table border="1"> <thead> <tr> <th></th> <th>茶渋</th> </tr> </thead> <tbody> <tr><td>O</td><td>54.1</td></tr> <tr><td>Ca</td><td>25.6</td></tr> <tr><td>Si</td><td>5.3</td></tr> <tr><td>Al</td><td>4.5</td></tr> <tr><td>Mg</td><td>2.9</td></tr> <tr><td>Cl</td><td>2.3</td></tr> <tr><td>S</td><td>1.6</td></tr> <tr><td>Mn</td><td>1.3</td></tr> <tr><td>P</td><td>1.2</td></tr> <tr><td>Na</td><td>0.7</td></tr> <tr><td>K</td><td>0.5</td></tr> </tbody> </table> <p>半定量値：wt%</p> <p>※主成分は有機物</p> | | 茶渋 | O | 54.1 | Ca | 25.6 | Si | 5.3 | Al | 4.5 | Mg | 2.9 | Cl | 2.3 | S | 1.6 | Mn | 1.3 | P | 1.2 | Na | 0.7 | K | 0.5 |  |
| | | 茶渋 | | | | | | | | | | | | | | | | | | | | | | | | |
| | O | 54.1 | | | | | | | | | | | | | | | | | | | | | | | | |
| | Ca | 25.6 | | | | | | | | | | | | | | | | | | | | | | | | |
| | Si | 5.3 | | | | | | | | | | | | | | | | | | | | | | | | |
| | Al | 4.5 | | | | | | | | | | | | | | | | | | | | | | | | |
| | Mg | 2.9 | | | | | | | | | | | | | | | | | | | | | | | | |
| | Cl | 2.3 | | | | | | | | | | | | | | | | | | | | | | | | |
| | S | 1.6 | | | | | | | | | | | | | | | | | | | | | | | | |
| | Mn | 1.3 | | | | | | | | | | | | | | | | | | | | | | | | |
| | P | 1.2 | | | | | | | | | | | | | | | | | | | | | | | | |
| | Na | 0.7 | | | | | | | | | | | | | | | | | | | | | | | | |
| K | 0.5 | | | | | | | | | | | | | | | | | | | | | | | | | |
| 備考 | | | | | | | | | | | | | | | | | | | | | | | | | | |

株式会社コーエキ — 総合分析と評価 —

長野県岡谷市田中町三丁目 3-24 TEL 0266-23-2155 FAX 0266-23-0733

e-mail : info@e-koeki.co.jp URL : <http://www.e-koeki.co.jp/index2.htm>